




REFERENCE MATERIAL PRODUCT INFORMATION SHEET

NMIA S045: Carboxy Finasteride

Report ID: S045.2024.01 (Ampouled 190711)

Chemical Formula: $C_{23}H_{34}N_2O_4$

Molecular Weight: 402.5g/mol

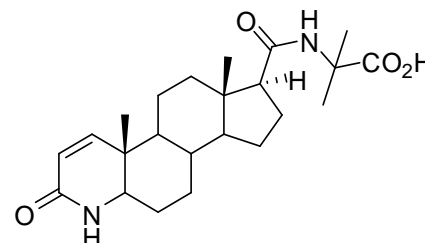

Property value

Batch No.	CAS No.	Estimated mass per ampoule
17-S-07	116285-37-1	820 ± 21 µg

The uncertainty has been calculated according to ISO Guide 35 and is stated at the 95% confidence limit ($k = 2$).

IUPAC name: N-[[[(4aR,4bS,6aS,7S,9aS,9bS,11aR)-4a,6a-Dimethyl-2-oxo-2,4a,4b,5,6,6a,7,8,9,9a,9b,10,11,11a-tetradecahydro-1H-indeno[5,4-f]quinolin-7-yl]carbonyl]-2-methylalanine.

Expiration of certification: The property values are valid till 16 October 2028, four years from the date of re-certification provided the **unopened** material is handled and stored in accordance with the recommendations below. The material as issued in the unopened container and stored as recommended below should be suitable for use beyond this date, subject to confirmation of batch stability from the issuing body. The expiry date/shelf life does not apply to ampoules that have been opened. In such cases it is recommended that the end-user conduct their own in-house stability trials.

Description: The compound is supplied as a dried aliquot in a sealed ampoule and is intended for a single use to prepare a standard solution containing S045. This material was prepared by sourced from an external supplier, and certified for identity and purity by NMIA.

Intended use: This reference material is recommended for qualitative analysis only and is not intended for use as a calibrator. The material does not have certified reference material status as metrological traceability of the stated purity value to the SI unit for mass (kg) has not been established.

Instructions for use: Open the ampoule and carefully rinse the interior at least three times with a suitable organic solvent (e.g. acetonitrile). This will transfer approximately 820 µg of anhydrous carboxy finasteride.

Recommended storage: When not in use, this material should be stored at or below 4 °C in a closed container in a dry, dark area.

Stability: In the absence of long term stability data the measurement uncertainty at the 95% coverage interval has been expanded to accommodate any potential change in the property value. The stability component has been estimated from stability trials conducted on similar materials by NMI Australia over the last ten years.

The long-term stability of the compound in solution has not been examined.

Homogeneity assessment: The homogeneity of the material was assessed using purity assay by HPLC with UV detection on seven randomly selected ampoules of the material. The material was judged to be sufficiently homogeneous at this level of sampling as the variation in analysis results between samples was not significantly different at a 95% confidence level from that observed on repeat analysis of the same sample.

Safety: Treat as hazardous substance. Use appropriate work practices when handling to avoid skin or eye contact, ingestion or inhalation of dust. Refer to the provided safety data sheet.

S. R. Davies

Dr Stephen R. Davies,
Team Leader,
Chemical Reference Materials, NMI.
21 October 2024

This report supersedes any issued prior to 21 October 2024.

NATA Accreditation No. 198 / Corporate Site No. 14214.

Legal notice: Terms and Conditions associated with the provision of this reference material can be found on the NMIA website.

Characterisation Report:

HPLC:	Instrument:	Waters Model 2695 / Thermofisher RS pump and auto sampler
	Column:	X-Bridge C-18, 5 µm (4.6 mm x 150 mm)
	Column oven:	40 °C
	Mobile Phase:	A = Milli Q water (0.05% Formic acid); B = Acetonitrile (0.05% Formic acid) 0-13 min 24% B; 13-20 min 24-80% B; 20-23 min 80%B; 23-24 min 80-24%B; 24-32 min 24% B
	Flow rate:	1.0 mL/min
	Detector:	Waters 2998 PDA / RS PDA, operating at 208 nm
	Relative peak area of the main component:	
	Initial analysis:	Mean = 89.6%, s = 0.3% (7 ampoules in duplicate, December 2019)
	Re-analysis:	Mean = 89.0%, s = 0.2% (5 ampoules in duplicate, November 2020)
	Re-analysis:	Mean = 89.5%, s = 0.1% (5 ampoules in duplicate, October 2021)
	Re-analysis:	Mean = 90.0%, s = 0.2% (5 ampoules in duplicate, September 2022)
	Re-analysis:	Mean = 91.0%, s = 0.02% (5 ampoules in duplicate, October 2024)

The following analytical data was obtained on the bulk material subsequently used in the preparation of the ampoules.

The identity was confirmed by a range of spectroscopic techniques, NMR, IR and MS. The certified purity value was obtained by mass balance from a combination of traditional analytical techniques, including HPLC with UV detection, thermogravimetric analysis, Karl Fischer analysis and ¹H NMR spectroscopy. The purity value is calculated as per Equation 1.

$$\text{Purity} = (100 \% - I_{\text{ORG}}) \times (100 \% - I_{\text{VOL}} - I_{\text{NVR}}) \quad \text{Equation 1}$$

I_{ORG} = Organic impurities of related structure, I_{VOL} = volatile impurities, I_{NVR} = non-volatile residue.

Supporting evidence is provided by qualitative elemental microanalysis.

HPLC:	Instrument:	Waters Model 1525 Binary pump, 717 plus autosampler
	Column:	X-Bridge C-18, 5 µm (4.6 mm x 150 mm)
	Column oven:	40 °C
	Mobile Phase:	A = MilliQ water; B = Acetonitrile 0-13 min 24% B; 13-20 min 24-80% B; 20-23 min 80%B; 23-24 min 80-24%B, 24-32 min 24% B Both aqueous and organic phases contained 0.05 % formic acid (v/v)
	Flow rate:	1.0 mL/min
	Detector:	Waters 2998 PDA operating at 208 nm
	Relative peak area of the main component:	
	Initial analysis:	Mean = 89.9%, s = 0.2% (7 sub samples in duplicate, January 2018)
	Re-analysis:	Mean = 91.0%, s = 0.4% (5 sub samples in duplicate, January 2019)
Karl Fischer analysis:		Moisture content 8.1% mass fraction (January 2018) Moisture content 9.0% mass fraction (January 2019)
Thermogravimetric analysis:		Volatile content 8.1% and non volatile residue 0.2 – 0.3 % mass fraction (April 2018)

Spectroscopic and other characterisation data

LC-MS:	<p>Instrument: Waters Acquity/Waters TQ Detector Column: X-Bridge C-18, 100 mm × 2.1 mm I.D. × 3.5 μm Column temp: 40 °C Solvent system: A = MilliQ water; B = Acetonitrile 0-13 min 24% B; 13-20 min 24-80% B; 20-23 min 80%B; 23-24 min 80-24% B, 24-32 min 24% B Both aqueous and organic phases contained 0.05 % formic acid (v/v) Flow rate: 0.2 mL/min Sample prep: 100 μg/g in mobile phase Injection volume: 10 μL Ionisation mode: Electrospray positive/negative ion Capillary voltage: 3.5 kV Cone voltage: 17 V Source temp: 120 °C Desolvation gas temp: 350 °C Cone gas flow rate: 1 L/hr Desolvation gas flow: 600 L/hr</p> <p>The retention time of carboxy finasteride is reported with the main peak in the mass spectrum. The latter is reported as a mass/charge ratio.</p> <p>4.3 min: 403.2 (M-H⁺) <i>m/z</i></p>
ESI-MS:	<p>Instrument: Micromass Quatro LC Micro Operation: Negative ion mode, direct infusion at 5.0 μL/min Ionisation: ESI spray voltage at 3.0 kV negative ion EM voltage: 650 V Cone voltage: 30 V Peak: 401.1 (M-H⁺) <i>m/z</i></p>
TLC:	<p>Conditions: Kieselgel 60F₂₅₄. Dichloromethane/methanol (85/15) Single spot observed, R_f = 0.31 (streaks) Visualisation with UV at 254 nm and permanganate</p>
IR:	<p>Instrument: Biorad FTS300MX FT-IR Range: 4000-400cm⁻¹, KBr powder Peaks: 3577, 3470, 3220, 3171, 3049, 2967, 2933, 2865, 2840, 1722, 1643, 1594, 1543, 1272, 1165, 813 cm⁻¹</p>
¹ H NMR:	<p>Instrument: Bruker Avance III-500 Field strength: 500 MHz Solvent: MeOH-<i>d</i>₄ (3.31 ppm) Spectral data: δ 0.71 (s, 3H), 0.96 (s, 3H), 1.05-1.10 (m, 2H), 1.17-1.22 (m, 1H), 1.26-1.64 (m, 12H), 1.66-1.74 (m, 3H), 1.77-1.84 (m, 2H), 1.97-2.01 (m, 1H), 2.15 (dd, 1H, <i>J</i>=10.8, 20.0 Hz), 2.26 (t, 1H, <i>J</i>=9.3 Hz), 3.33-3.35 (m, 1H), 5.76 (d, 1H, <i>J</i>=9.9 Hz), 6.97 (d, 1H, <i>J</i>=9.9 Hz), 7.92 (s, 1H) ppm Both dichloromethane and acetone were measured at < 0.1% mass fraction, respectively.</p>
¹³ C NMR:	<p>Instrument: Bruker DMX600 Field strength: 151 MHz Solvent: MeOH-<i>d</i>₄ (49.0 ppm) Spectral data: δ 12.2, 13.9, 22.2, 24.4, 25.3, 25.4, 25.7, 26.3, 30.7, 36.6, 38.8, 40.4, 45.9, 56.86, 56.94, 57.0, 60.9, 123.1, 153.6, 168.8, 174.5, 178.2 ppm</p>
Microanalysis:	<p>Found: C = 62.5%; H = 9.0%; N = 6.4% (January 2018) Calculated: C = 63.0%; H = 8.7%; N = 6.4% (Calculated for C₂₃H₃₄N₂O₄·2H₂O)</p>